


SI Digital PC+ 'S'

New Generation Table Top Sonic Irrigator

The SI Digital PC+ 'S' gives consistent and unequalled cleaning results every time. The SI Digital PC+ 'S' offers Sonic Irrigation with CPE® (Cannulated Pulse Enhancement) with the addition of monitoring features to record key performance data for instrument traceability.

Key benefits & features

- Sonic Irrigation™
- 28L tank total capacity
- Process up to 12 cannulated instruments
- New touch screen & animation technology
- USB port connectivity – Trackability & Backup
- E-Zee load tray system compatibility
- Process Interlocked lid
- CPE® (Cannulated Pulse Enhancement)
- Load performance optimisation
- Pre settable process time / temperature
- HTM 01-01 compliant
- UL & CE compliant

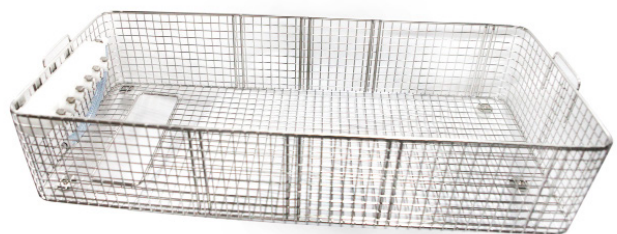


Unique elements




Fast Cycle Time



Touch Screen / Animated Graphic Interface

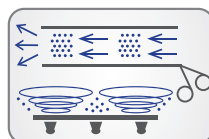

Sonic Irrigation

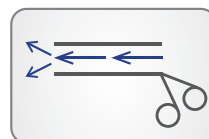

Lumen Flush



RS232 Output Data connection for Printer

SI Digital PC+ 'S'

New Generation Table Top Sonic Irrigator

Medisafe

The SI Digital PC+ 'S' system monitors the critical cleaning process parameters such as time, temperature and ultrasonic power in accordance with current guidance and international standards. It can be recorded permanently through a discreetly placed USB port which intelligently detects whether you wish to keep the information on a database device for traceability and legislative issues.

Machine includes

- EZ Load Basket
- SI Digital Starter Plumbing Kit
- Instruction for use

Optional Accessories

- M10184 Dosing Pump

Associated Detergents


	4L	5L	10L
3E-Zyme	MED8035	MED8010	MED8037
3E-Zyme Scope Plus	M20029	MED8015	M20200
4-Zyme		MED8020	MED8021
EDA+		MED7700	MED7701
Medi-pH- Safe		M20011	M20008

Machine specification

Use: For internal and external cleaning of fine surgical instruments and other lumen instruments

Dimensions (external):

Height: 345mm (13.5")
Width: 715mm + drain 40mm (28" + 1.5")
Depth: 445mm (17.5")

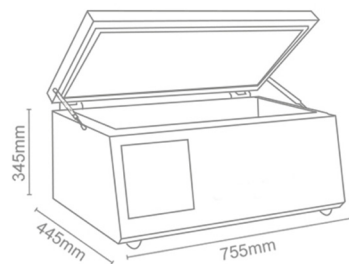

Electrical Power:

1200kW (in use)

Tank Total Volume:

28L

Water Supply:

Manual filling (not automatic)

Environmental Conditions:

For indoor use at altitudes up to 2000m. Operating temperature: 10-40 °C. 80% for temperature up to 31 °C decreasing to 50% at 40°C. Designed for pollution degree 2, i.e. not totally sealed.

Unit weight:

30kg (66lbs) empty

Capacity:

Up to 12 cannulated instruments,
General Surgical Instruments

Electrical Supply:

100-127 Vac 50Hz or 60Hz
220-240 Vac 50Hz or 60Hz

Construction:

High grade 316 (Tank), 304 (general construction) corrosion resistant stainless steel with brushed finish.

Order details

MED521035 - 220-240 Vac-50Hz

MED521036 - 220-240 Vac-60Hz

MED521020 - 100-127 Vac-50/60Hz